
Product Sheet



Catalog Number

9502C-S34-T-M

Caps, Snap Top, with 1.00mm thick, White Silicone Rubber / Red PTFE Pre-Slit (Y Shape) Septa, fitted in smooth, soft polyethylene, Natural Transparent Caps. For use with 11mm, 12x32mm Snap Top autosampler

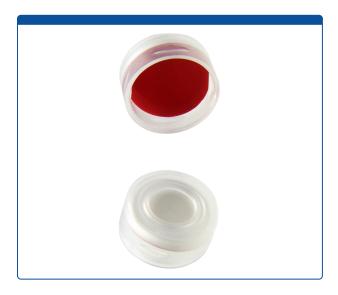

Additional Info

Soft polyethylene caps are more pliable than other Snap Caps and are easier to place on the vial. Popular choice when using a large number of caps at one time. These caps may not seal as well as other Snap Caps and should be used with sample solvents with low volatility.

Y Shape Pre-Slit provides a greater target area for needles and are easier to penetrate. The 1.0mm Septa is thicker than most and offers improved sealing compared with other similar Septa.

Pre-Slit Septa are normally recommended for Waters Alliance and Acquity instruments.

Click **HERE** for smaller pack quantities.

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451 tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

\$1,271.44

Printed On 2024/05/20